



MEDICAL AND REHABILITATION  
EQUIPMENT

[www.innow.eu](http://www.innow.eu)



### **EXPERIENCE**

on the market since 1984

### **QUALITY**

tens of thousands of users

### **SENSITIVITY**

solutions to every need

### **SATISFACTION**

satisfaction is a priority



# **TREATMENT CHAIRS AND TABLES**

## **MANUFACTURER'S CATALOGUE**

# TREATMENT CHAIRS AND TABLES



Treatment chairs and tables are widely used in many areas of medicine: starting from ophthalmology and laryngology and ending up with surgery and healthcare cosmetology.



*In the photo:  
FZ-02 GINN AT gynaecological chair  
with optional accessories: anatomical  
upholstery, medical paper roll holder.*

Many years of experience, professionalism, CE marking, high quality, and continuous development of products allow us to adapt these products to the needs of patients and medical personnel. A wide range of upholstery colours ensures the adaptation of treatment chairs and tables to individual user preferences.

The manufacturing materials are durable and easy to clean.



*In the photo:  
S412E [KARDIO]  
treatment (cardiological) table  
with additional  
accessories: side barrier,  
medical paper roll holder.*



*In the photo:  
ZW100 treatment table with optional accessories:  
Trendelenburg and anti-Trendelenburg position,  
side rails for installing additional accessories,  
medical paper roll holder.*

All products come with washable upholstery whose colour can be selected from several available in the offer.





seat/lying surface  
height



seat/lying surface  
width



seat/lying surface  
depth



maximum  
load



seat tilt angle



backrest  
tilt angle

## FOL-01

### Ophthalmology and ENT chair

A simple treatment chair designed for ophthalmological and laryngological procedures. It is controlled hydraulically. The seat is fastened to an aesthetically-pleasing chrome-plated base. It comes with a comfortable headrest for head stabilisation with height, extension, and tilt angle adjustment. The chair also has backrest tilt angle adjustment, increasing user comfort even more.

#### Optional accessories:

1. multi-position polyurethane armrest



490–620



470



430



125



20



## FOL-02

### Ophthalmology and ENT chair

A multifunctional treatment chair with a 3-section lying surface located on a chrome-plated and cone-shaped base. Height adjustment is possible thanks to a hydraulic cylinder. The backrest and leg segment are controlled by gas springs, which make it possible for the chair to move to sitting, lying, and intermediate positions. The chair comes with an upholstered headrest adjustable to many positions. It is used for ophthalmological and laryngological procedures.

#### Optional accessories:

1. standard upholstered armrest
2. multi-position polyurethane armrest
3. merged upholstery (two-coloured)
4. option to lower the backrest and lift the footrest simultaneously and vice versa



550–700



500



500



150






FZ01 EU [3CP]

**EUREKA treatment chair with seat tilt angle adjustment**

A modern treatment chair with a 3-section lying surface: backrest, seat, leg part. The position of the backrest and the leg part is controlled by a gas spring. The chair can go to sitting, lying, and intermediate positions. Seat tilt angle adjustment makes it possible for the chair to change to cradle, catapult, Trendelenburg, and anti-Trendelenburg positions.

**Optional accessories:**

1. side holder
2. two-section leg part
3. extendable leg part
4. side rails for installing additional accessories
5. medical paper roll holder
6. tiltable headrest with height adjustment
7. pillow-type headrest
8. headrest with height adjustment
9. footrest
10. standard upholstered armrest
11. standard polyurethane armrest
12. multi-position polyurethane armrest
13. anatomical upholstery
14. drip holder

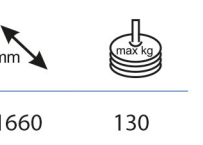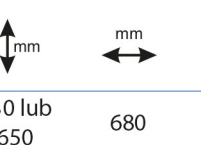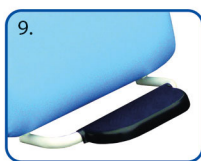

mm  
530 lub  
650

mm  
680

mm  
1660

max kg  
130






**EUREKA treatment chairs**

Modern treatment chairs with a 3-section lying surface: backrest, seat, leg part. The position of the backrest is controlled by a gas spring. The chair can go to sitting, lying, and intermediate positions.

**Optional accessories:**



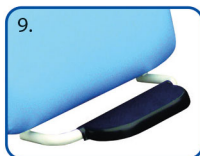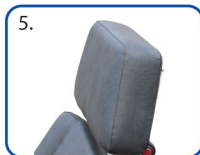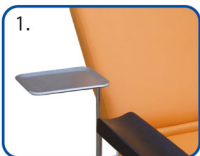

1. side holder
2. two-section leg part
3. extendable leg part
4. side rails for installing additional accessories
5. medical paper roll holder
6. tiltable headrest with height adjustment
7. pillow-type headrest
8. headrest with height adjustment
9. footrest
10. standard upholstered armrest
11. standard polyurethane armrest
12. multi-position polyurethane armrest
13. anatomical upholstery
14. drip holder



FZ01 EU [3C]



FZ01 EU [3CW]  
with height  
adjustment

|                      |  mm |  mm |  mm |  max kg |
|----------------------|--|--|--|--|
| <b>FZ01 EU [3C]</b>  | 530 lub<br>650   | 680  | 1660   | 130  |
| <b>FZ01 EU [3CW]</b> | 530–830  | 680  | 1660   | 130  |

## FZ01 EU [2C], FZ01 EU [2CW]

## EUREKA treatment chairs

Ergonomic treatment chairs with a 2-section lying surface and robust design. The backrest is controlled manually with gas springs. Sitting, intermediate, and lying positions are available.

## Optional accessories:





1. side holder
2. side rails for installing additional accessories
3. medical paper roll holder
4. tiltable headrest with height adjustment
5. pillow-type headrest
6. headrest with height adjustment
7. standard upholstered armrest
8. standard polyurethane armrest
9. multi-position polyurethane armrest
10. anatomical upholstery
11. drip holder



FZ01 EU [2C]



FZ01 EU [2CW] with height adjustment

|                      |  mm |  mm |  mm |  max kg |
|----------------------|--|--|--|--|
| <b>FZ01 EU [2C]</b>  | 530  | 680  | 1300   | 130  |
| <b>FZ01 EU [2CW]</b> | 530–830  | 680  | 1300   | 130  |

### Treatment chair

A multifunctional treatment chair designed for procedures, collections, and dialysis, as well as for transport and care of the patient after such procedures. It is used in laryngology or cardiology, as well as during chemotherapy. It comes with a 3-section upholstered lying surface, which can be adjusted to sitting, lying, and intermediate positions. There are four braking casters for easier chair transport. Lifting mechanisms are protected with an aesthetically-pleasing plastic cover.




FZ02 - basic chair version

FZ02 [W] - version with height adjustment

FZ02 [P] - version with seat tilt angle adjustment

### Optional accessories:

1. side holder
2. two-section leg part
3. extendable leg part
4. side rails for installing additional accessories
5. medical paper roll holder
6. tiltable headrest with height adjustment
7. pillow-type headrest
8. headrest with height adjustment
9. footrest
10. standard upholstered armrest
11. standard polyurethane armrest
12. multi-position polyurethane armrest
13. anatomical upholstery
14. drip holder

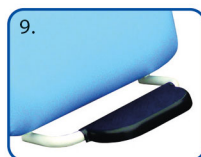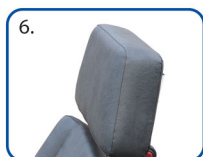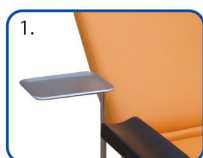

**FZ02, FZ02 [P]**

600

600

480

130

**FZ02 [W]**

600-880

600

480

130



## F301 EVO, F301H EVO

**TRAPER II rehabilitation, care, and geriatric chair for:**

- the elderly,
  - the disabled,
  - people who need constant care,
- for transporting people both inside and outside buildings.

F301 EVO and F301H EVO models are equipped with large front wheels for outdoor transport – on an uneven surface as well.



F301H EVO  
with height  
adjustment

**Standard equipment:**

- **synchronised adjustment** – fast and simultaneous backrest lowering and leg part lifting to obtain the lying position (and vice versa – the sitting position),
- **hydraulic seat tilt angle adjustment** (F301H EVO and F301KL H EVO models) – makes it possible for the chair to go to the cradle position that protects the patients from slipping off the chair (max. angle 24°) and to the catapult position that makes it easier for the patient to get off the chair (max. angle 5°),
- **footrest** – when the weight of the body is on it, it becomes level and lowers to the ground, making it easier for the patient to get up from the chair; once the weight is off the footrest, it moves to the original position (CROSS version with a removable footrest for outdoor use),
- **upholstered and removable armrests** with height adjustment,
- **modern design** that makes it easy and convenient for medical personnel to give care treatment,
- **easy adjustment** – the position of the patient's body can be changed frequently, lowering the risk of bedsores,
- **highest-quality materials** – the chair is comfortable for its users,
- **many upholstery colours to choose from** – ideal adaptation of the chair to individual preferences.



F301 EVO  
CROSS



F301 EVO

**Examples of F301H EVO table positions**

## F301KL H EVO

### TRAPER II postoperative, care, geriatric chair for:

- patients after procedures and surgeries,
- the elderly,
- the disabled,
- people who need constant care, to be used inside buildings.

F301KL EVO and F301KL H EVO models have swivel casters for easy movement in narrow passages and tight spaces.

F301KL H EVO,  
with tilt angle  
adjustment



Optional accessories for all models (F301 EVO, F301H EVO, F301KL EVO, F301KL H EVO):

1. removable wooden table
2. seat belt
3. five-point seat belt
4. pillow-type headrest
5. profiled pillow-type headrest
6. stabilising headrest with height adjustment
7. headrest with height adjustment
8. anatomical upholstery
9. drip holder
10. walking cane holder
11. cup holder
12. removable CROSS footrest



F301KL EVO



500



500



495



150

## FZ-02 GINN

**Gynaecological chair**

A functional treatment chair for gynaecological and urological examinations and procedures. It has a 2-section lying surface with an opening. Backrest tilt angle adjustment with a gas spring. Electrical lying surface height adjustment with a foot remote control (hand remote control – option). The set includes 2 fixed armrests, 2 multi-position polyurethane thigh supports, a chassis cover, and a gynaecological bowl.

**Optional accessories:**

1. 2 rollers for moving the chair
2. side rails for installing additional accessories on the seat segment
3. hand remote control
4. medical paper roll holder
5. headrest adjustable in the range of 0-45°
6. multi-position polyurethane armrest
7. anatomical upholstery



FZ-02 GINN with optional accessories: side rails, hand remote control.



530-960



690



1425



130



### Gynaecological chair with Trendelenburg and anti-Trendelenburg positions

A functional treatment chair for gynaecological and urological examinations and procedures. It has a 2-section lying surface with an opening. Backrest tilt angle adjustment with a gas spring. Electrical lying surface height adjustment with a foot remote control (hand remote control – option), The set includes 2 fixed armrests, 2 multi-position polyurethane thigh supports, a chassis cover, and a gynaecological bowl. The FZ-02 GINN AT model can also be moved to Trendelenburg (10°) and anti-Trendelenburg (15°) positions with a foot remote control.

#### Optional accessories:

1. 2 rollers for moving the chair
2. side rails for installing accessories on the seat segment
3. hand remote control
4. medical paper roll holder
5. headrest adjustable in the range of 0-45°
6. multi-position polyurethane armrest
7. anatomical upholstery



FZ-02 GINN AT with optional accessories: anatomical upholstery, medical paper roll holder.



575–1000



640



1425



130

## FG100

**Gynaecological table**

A gynaecological table with a robust steel base, equipped with rollers for easy transport. Height adjustment using a hydraulic cylinder and two levers located on both table sides. Backrest adjustment using a gas spring. The set includes a gynaecological bowl.

**Optional accessories:**

1. stainless steel base
2. side rails for installing additional accessories
3. medical paper roll holder
4. pillow-type headrest
5. multi-position polyurethane armrest
6. multi-position polyurethane thigh support
7. Trendelenburg and anti-Trendelenburg positions



FG100 with optional accessories: side rails for installing additional accessories, medical paper roll holder, multi-position polyurethane thigh support, Trendelenburg and anti-Trendelenburg positions



FG100 with optional accessories:  
multi-position polyurethane thigh supports  
side rails for installing accessories, medical paper roll holder.



700–1100



600



1500



150





### Treatment table

A simple treatment table with height adjustment by steps. A robust powder-coated base, additionally equipped with a levelling foot. A two-section lying surface with an opening and a gynaecological bowl included. The table is designed for performing procedures and examinations in various areas of medicine.

#### Optional accessories:

1. gas spring lifted backrest
2. lying surface with a headrest – extends the table to 175 cm
3. side rails for installing additional accessories
4. medical paper roll holder
5. pillow-type headrest
6. multi-position polyurethane armrest
7. multi-position polyurethane thigh support
8. frame for holding anaesthesia screen
9. split headrest – option available when the backrest is adjusted mechanically

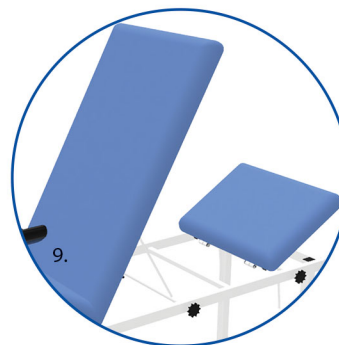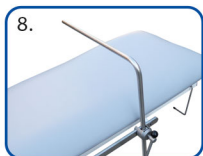

800 – 850 –  
900



600



1500



150



## FZ-02 FINN, FZ-02 FINN AT



## Treatment table





A multifunctional treatment table used in surgery, urology, and gynaecology. The basic version comes with a 4-section lying surface. Electrically adjustable height with a hand remote control (foot remote control – option), backrest adjustment with a gas spring. The leg part (split into two parts) is also adjusted with a gas spring. The set also includes a lower frame cover and side rails on the backrest, seat, and leg part segment for installing additional accessories.

## Optional accessories:

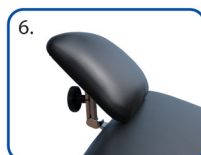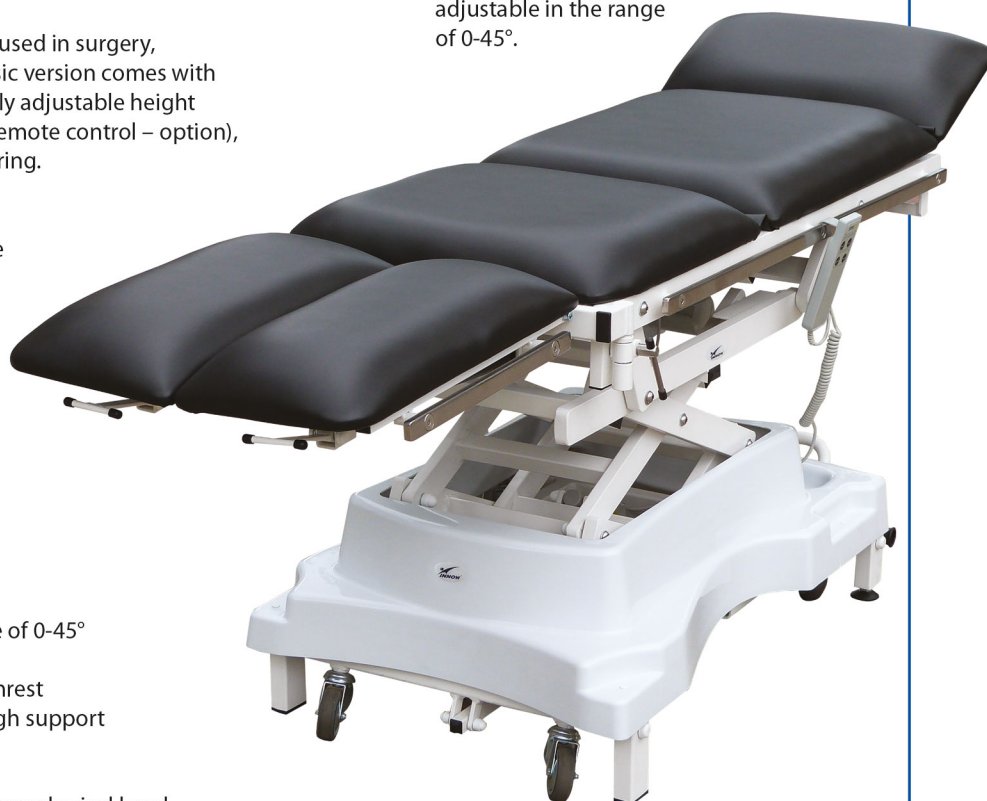1. 2 rollers for moving the table
2. lifting section cover
3. foot remote control
4. medical paper roll holder
5. headrest adjustable in the range of 0-45°
6. multi-position headrest
7. multi-position polyurethane armrest
8. multi-position polyurethane thigh support
9. stainless steel frame
10. seamless upholstery
11. opening in upholstery and gynaecological bowl

FZ-02 FINN AT with Trendelenburg and anti-Trendelenburg positions and optional accessories: headrest adjustable in the range of 0-45°.



|                       |  mm |  mm |  mm |  max kg |
|-----------------------|--|--|--|--|
| <b>FZ-02 FINN</b>     | 605–875  | 550  | 1975   | 130  |
| <b>FZ-02 FINN AT*</b> | 700–970  | 550  | 1975   | 130  |



FZ-02 FINN

\*Examples of FZ-02 FINN AT table positions





### Treatment table

A well-equipped and multifunctional treatment table with Trendelenburg and anti-Trendelenburg positions as well as washable upholstery for medical applications. It comes with: a lying surface with a separate backrest, seat, leg part; a multifunctional headrest for easy patient head adjustment during procedures, a hand remote control for height and Trendelenburg and anti-Trendelenburg position adjustment.

The table also features: backrest adjustment with a gas spring, lower frame cover, side rails on the lying surface and leg part segment for installing additional accessories, lifting section cover, upper frame with a cover plate, and a medical paper roll holder.

### Optional accessories:

1. 2 rollers for moving the table
2. foot remote control
3. multi-position polyurethane armrest
4. multi-position polyurethane thigh support
5. stainless steel frame
6. seamless upholstery
7. opening in upholstery and gynaecological bowl



700–970



550



1975



130





## ZW100

## Treatment table

A simple treatment table with a scissor-shaped frame. It comes with a 5-section lying surface with an opening and a gynaecological bowl. It has 2 rollers for easy transport. Height adjustment with a hydraulic cylinder and two side levers. Backrest tilt angle adjustment with a gas spring. The table can be quickly converted into a gynaecological chair by lifting the backrest, sliding the leg part to the side, and adding thigh supports.

## Optional accessories:

1. stainless steel base
2. side rails for installing additional accessories
3. side rails for installing additional accessories on the leg support segment
4. medical paper roll holder
5. multi-position polyurethane armrest
6. multi-position polyurethane thigh support
7. Trendelenburg (10°) and anti-Trendelenburg (10°) positions
8. frame for holding anaesthesia screen
9. leg support segment adjustment
10. seamless upholstery



630–1050



600



2000



150

750–1170

600

2000

150



**Electric cardiological table**

It is designed for cardiological procedures and procedures performed with the patient lying on the side. It has a special opening in the lying surface and a side drawer. The height-adjustable headrest for procedures performed with the patient lying on the side ensures that the patient can be precisely positioned as required.

**Optional accessories:**

1. side barrier
2. side rails for installing additional accessories
3. medical paper roll holder
4. hydraulic height adjustment (S412H [KARDIO])
5. emergency switch



S412E [KARDIO] with optional accessories:  
side barrier, medical paper roll holder.



500–900



700



2000



200

# TABLE ACCESSORIES

## Stabilising clamps:



**OMINN3900**  
Standard clamp



**FINN176**  
Multi-position clamp



**OMINN4050**  
Multi-position universal  
clamp



**CHROFINN-STZ01**  
Economic normal clamp  
(chrome-plated)

## Armrests and thigh supports – with bell joints:



**FINN055**



**OMINN2200**



**OMINN1200**



**FINN075**

### Auxiliary equipment:



**FINN 015**  
Frame for holding  
anaesthesia  
screen installed on  
side rail



**OMINN3650**  
Arm support



**OMINN1100**  
Drip holder  
installed on  
side rail



**FINN325**  
Side barriers

### Armrests and thigh supports – standard:



**PODINN\_arm**  
Multi-position  
polyurethane  
arm rest



**PODINN\_thigh**  
Multi-position  
polyurethane thigh  
support

### Stabilising belts:



**OMINN1550**



**OMINN1500**



**OMINN1520**

### Gel cushions:



**OMINN4100**



**OMINN4200**





#### **SALES OFFICE**

**40-763 Katowice, ul. Ustrońska 8  
POLAND  
www.innow.eu, on-line 24h**

#### **SALES SHOWROOM**

**40-763 Katowice, ul. Ustrońska 8  
open: Mon. – Fri. 8.00 am – 4.00 pm  
Medical furniture is on display permanently and  
can be purchased in the showroom directly.**

#### **CONTACT**

**tel. +48 32 252 67 56  
salon@innow.pl**